

Concise Communication

In vitro assessment of effect of initial specimen diversion device on detection of central venous catheter contamination or colonization

Mark E. Rupp MD¹ , Paul D. Fey PhD², Elizabeth Lyden MS³  and Luke Handke PhD²

¹Division of Infectious Diseases, University of Nebraska Medical Center, Omaha, NE, USA, ²Department of Pathology, Microbiology, and Immunology, University of Nebraska Medical Center, Omaha, NE, USA and ³Department of Biostatistics, University of Nebraska Medical Center, Omaha, NE, USA

Abstract

The role of initial specimen diversion devices (ISDDs) in preventing contamination of central venous catheter (CVC) blood cultures is undefined. A model to simulate CVC colonization and contamination compared standard cultures with ISDD technique. ISDD detected 100% of colonized CVCs while decreasing false-positive cultures from 36% to 16%.

(Received 18 September 2024; accepted 20 November 2024; electronically published 19 December 2024)

Introduction

Blood cultures are commonly used and critically important.¹ 1%–5% of blood cultures are contaminated by skin-residing commensal organisms, resulting in adverse outcomes.^{1–3}

Catheter-related bloodstream infections (CRBSIs) are common and result in morbidity and excess cost.^{1–4} Diagnosis of CRBSI often requires cultures from the vascular catheter and peripheral blood.⁵ Cultures obtained from vascular catheters are more likely than peripheral cultures to be contaminated,⁶ often leading to unnecessary catheter removal.

Initial specimen diversion devices (ISDDs) cost-efficiently decrease blood culture contamination.⁷ It is unknown whether ISDD technology is helpful in decreasing contamination of central venous catheter (CVC)-drawn cultures. We developed in vitro models of CVC contamination and colonization to characterize the utility of the ISDD in detecting CVC colonization while minimizing culture contamination.

Methods: (Also see Supplemental Materials)

CVC colonization

CVC (Arrow triple-lumen, Teleflex) colonization was simulated by filling the lumen with Trypticase Soy Broth (TSB) containing 10 colony-forming units (CFUs) of *S. epidermidis* 1457/pCM29 (*ica*-positive, biofilm-producing, expressing green fluorescent protein (GFP)). The CVCs, maintained in sterile catheter shields (Cath-Gard, Teleflex), were incubated overnight at 37°C. 50 CVCs were sampled by standard or ISDD method.

Corresponding author: Mark E. Rupp; Email: merupp@unmc.edu

Cite this article: Rupp ME, Fey PD, Lyden E, Handke L. In vitro assessment of effect of initial specimen diversion device on detection of central venous catheter contamination or colonization. *Infect Control Hosp Epidemiol* 2025. 46: 197–200, doi: [10.1017/ice.2024.220](https://doi.org/10.1017/ice.2024.220)

© The Author(s), 2024. Published by Cambridge University Press on behalf of The Society for Healthcare Epidemiology of America. This is an Open Access article, distributed under the terms of the Creative Commons Attribution licence (<https://creativecommons.org/licenses/by/4.0/>), which permits unrestricted re-use, distribution and reproduction, provided the original article is properly cited.

Catheter contamination

Catheter connector (MaxZero™, Becton Dickinson) contamination was simulated by inoculation of 10 CFUs of *S. epidermidis* in 50 µL of TSB onto the diaphragm of the CVC connector and allowing it to air dry. Fifty CVCs were sampled by standard or ISDD method.

Blood culture

Blood cultures were simulated by drawing sterile phosphate-buffered saline (PBS) through the CVC. In the standard method, the ISDD (SteriPath Blood Collection System, Magnolia Medical Technologies) was pre-engaged (diversion chamber closed) and then connected to the CVC. For the ISDD method, the ISDD was engaged to divert the initial 1.5–2 mL of PBS. For both methods, a Vacutainer (Becton Dickinson) was used to collect 1 mL of PBS for quantitative culture. Two 10 mL samples were then collected into BACTEC Plus Aerobic/F blood bottles (Becton Dickinson). The bacterial titer in CFU/mL was determined. Bottles were incubated in a BACTEC FX instrument (Becton Dickinson) and monitored for 5 days. Time-to-positivity (TTP) was recorded.

Dose-ranging preliminary study

To establish inoculum levels to simulate CVC connector contamination and catheter colonization, a dose-ranging study was conducted using 3 CVCs per inoculum group: 10¹ CFU–10³ CFU in colonization model; 10¹ CFU–10⁴ CFU in contamination model.

Scanning electron microscopy (SEM)

(See supplemental materials).

Statistical analysis

Fisher's exact test compared the proportion of positive cultures between the ISDD and standard methods. Wilcoxon rank sum test compared the median CFU/mL between the blood culture methods.



Table 1. Dose-ranging study. Quantitative culture results and blood culture time-to-positivity for colonization model and contamination model

Inoculum 10 ¹ CFU	Blood culture method	
	ISDD	Standard
Colonization model		
Quantitative blood culture (mean CFU/mL + SD)	1.4 × 10 ⁵ ± 3.4 × 10 ⁴	1.8 × 10 ⁷ ± 6.1 × 10 ⁶
1 st bottle TTP (hours) (mean + SD)	10.97 ± 0.93	7.75 ± 0.33
2 nd bottle TTP (hours) (mean + SD)	12.17 ± 0.39	10.2 ± 0.45
Inoculum 10² CFU		
Quantitative blood culture (mean CFU/mL + SD)	8.6 × 10 ⁶ ± 3.7 × 10 ⁶	9.5 × 10 ⁷ ± 3.7 × 10 ⁷
1 st bottle TTP (hours) (mean + SD)	7.7 ± 1.0	6.1 ± 0.2
2 nd bottle TTP (hours) (mean + SD)	10.0 ± 1.6	8.24 ± 0.7
Inoculum 10³ CFU		
Quantitative blood culture (mean CFU/mL + SD)	4.6 × 10 ⁷ ± 1.2 × 10 ⁷	3.0 × 10 ⁸ ± 3.0 × 10 ⁷
1 st bottle TTP (hours) (mean + SD)	5.06 ± 0.0	5.0 ± 0.1
2 nd bottle TTP (hours) (mean + SD)	6.17 ± 0.6	6.0 ± 0.8
Contamination model		
Quantitative blood culture (mean CFU/mL + SD)	0 ± 0	10 ± 10
1 st bottle TTP (hours) (mean + SD)	22.4 (1 positive)	22.1 ± 0.7 (2 positive)
2 nd bottle TTP (hours) (mean + SD)	0 (0 positive)	0 (0 positive)
Inoculum 10² CFU		
Quantitative blood culture (mean CFU/mL + SD)	10 ± 10	16.7 ± 5.8
1 st bottle TTP (hours) (mean + SD)	21.1 ± 0.5 (3 positive)	19.9 ± 0.3 (3 positive)
2 nd bottle TTP (hours) (mean + SD)	0 (0 positive)	0 (0 positive)
Inoculum 10³ CFU		
Quantitative blood culture (mean CFU/mL + SD)	50 ± 86	180 ± 112
1 st bottle TTP (hours) (mean + SD)	18.7 ± 0.3 (3 positive)	18.7 ± 0.6 (3 positive)
2 nd bottle TTP (hours) (mean + SD)	21.5 ± 0.7 (3 positive)	20.9 ± 0.5 (3 positive)
Inoculum 10⁴ CFU		
Quantitative blood culture (mean CFU/mL + SD)	343 ± 421	7.1 × 10 ³ ± 9.6 × 10 ³
1 st bottle TTP (hours) (mean + SD)	17.9 ± 0.7 (3 positive)	16.4 ± 0.4 (3 positive)
2 nd bottle TTP (hours) (mean + SD)	19.6 ± 0.1 (3 positive)	18.7 ± 0.1 (3 positive)

Note. CFU, colony-forming units; ISDD, initial specimen diversion device; TTP, time-to-positivity.

Results

Dose-ranging study

Colonization: All inoculum levels (10¹–10³ CFU) resulted in 100% blood culture positivity. Table 1 details quantitative culture and TTP data. For the full study, the lowest inoculum (10¹ CFU) that reliably resulted in colonized CVCs was used.

Contamination: At the 10¹ level, a differentiation was evident between ISDD and standard techniques. At larger inoculum (10²–10⁴ CFU), both the ISDD and standard cultures were positive. For the full study, the 10¹ CFU level was chosen.

Full study

Catheter colonization: All simulated blood cultures yielded bacterial growth in the standard and ISDD arms. There was a lower number of bacteria in the ISDD cultures versus standard cultures (1.17 × 10⁵ CFU/mL ± 1.85 × 10⁵ (SD) vs 3.96 × 10⁶ CFU/mL ± 4.97 × 10⁶,

respectively) ($P < 0.0001$). There was a significant difference in TTP for the initial bottle (12.51 h ± 1.61 h vs 10.02 h ± 1.03 h) ($P < 0.0001$) and the second bottle (14.14 h ± 1.81 h vs. 12.92 h ± 1.18 h) ($P = 0.017$) for ISDD and standard cultures, respectively.

Connector contamination: 4 of 25 (16%) ISDD cultures and 9 of 25 (36%) standard cultures yielded bacteria ($P = 0.196$). 0 of 25 quantitative ISDD cultures recovered bacteria, and 2 of 25 standard cultures yielded growth. There was not a significant difference in TTP between ISDD and standard cultures. First bottle positivity was 23.76 h ± 1.83 h versus 22.94 h ± 1.2 h ($P = 0.247$) for the ISDD and standard cultures, respectively. Second bottle positivity was 28.4 h (only 1 positive bottle) versus 22.56 h ± 1.26 h ($P = 0.54$) for the ISDD and standard cultures, respectively.

Inadvertent contamination was excluded with all arising colonies demonstrating GFP production.

Four colonized CVCs and two contaminated connectors were examined by scanning electron microscopy (SEM). All colonized catheters exhibited widespread adherent staphylococci (Figure 1A),

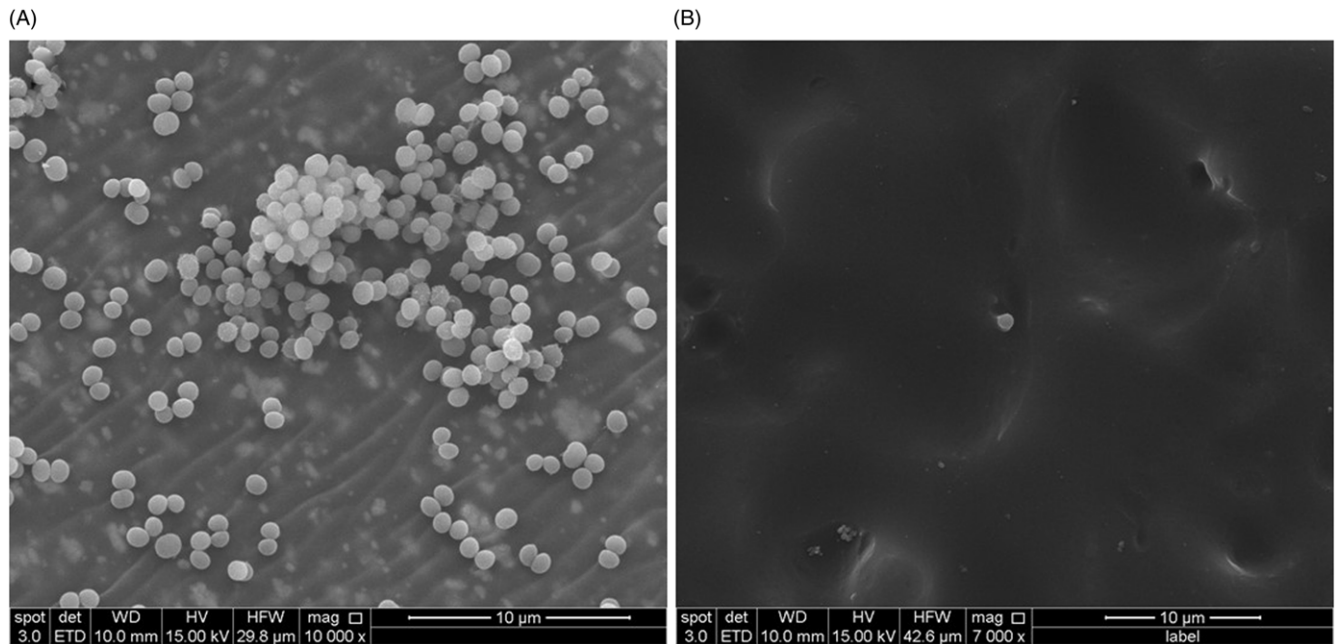


Figure 1. (A) Representative scanning electron micrograph of colonized catheter demonstrating numerous intraluminal adherent staphylococci. (B) Representative scanning electron micrograph of contaminated diaphragmatic surface of CVC connector valve demonstrating rare adherent staphylococci.

while only rare bacteria were visualized on the contaminated connectors (Figure 1B).

Discussion

ISDD prevents contamination of peripheral blood cultures.⁷ CVC blood cultures are more likely to be contaminated than peripheral cultures.⁶

Whether ISDD limits contamination of CVC cultures is unknown. Because the ISDD diverts the initial portion of the specimen (excluding blood with the heaviest burden of bacteria), the ISDD could result in exclusion of colonized catheters. Ideally, use of the ISDD would not decrease sensitivity while improving specificity (prevention of contamination due to bacteria on the connector). This *in vitro* study, modeling CVC intraluminal colonization and connector contamination, suggests possible clinical utility for ISDD CVC cultures.

In the colonization study, standard cultures and ISDD cultures detected 100% of colonized CVCs. Because the ISDD excludes the initial more heavily colonized portion, there was a longer TTP for ISDD cultures compared to standard cultures (2.49 hours longer incubation for first bottle, 1.22 hours longer for second bottle). This is supported by the greater number of bacteria present in quantitative cultures (1.17×10^5 CFU/mL (ISDD) vs 3.96×10^6 CFU/mL (standard)). The longer TTP noted with ISDD cultures would generally not be regarded as clinically significant. Thus, the ISDD method detects colonized CVCs without loss of sensitivity.

The utility of the ISDD technology for CVC cultures comes with improved specificity. The contamination study, which modeled low-level contamination of CVC connectors, demonstrated a decrease in blood culture contamination from 36% to 16% (56% relative risk reduction). This observation did not reach the level of statistical significance ($P=0.196$) most likely due to study size/statistical power. The low inoculum level for the contamination model is supported by clinical studies examining connector contamination in which most contaminated connectors exhibited

only a few CFU.¹⁰ At our center, despite robust infection prevention and stewardship programs, blood culture contamination results in an extra day of hospitalization and antibiotic treatment.² Additional adverse effects include unnecessary laboratory studies, unneeded antibiotics, diagnostic confusion, and unnecessary removal of CVCs.^{3,8,9}

The CVC colonization model resulted in a heavily colonized intraluminal surface from a small inoculum (10 CFU) after a short incubation time (24 hours) in the presence of a rich nutritional environment. This is analogous to a CVC used to instill total parenteral nutrition, lipids, or blood products. The contamination model mimics the small inoculum associated with contamination of the catheter hub or connector and less frequently resulted in positive cultures (36%) and rarely was detected by quantitative culture (8%). The dose-ranging study indicated that at higher levels of connector contamination, the discriminatory effect of the ISDD was lost. However, we believe low-level inoculum is a more accurate reflection of clinical conditions.¹⁰

In this CVC contamination/colonization model, the ISDD resulted in 100% sensitivity and improved specificity and justifies a prospective clinical trial.

Supplementary material. The supplementary material for this article can be found at <https://doi.org/10.1017/ice.2024.220>.

Data availability statement. All raw data are available upon request.

Acknowledgments. None.

Author contributions. MER, LH, and PDF contributed to the conception and design of the study. MER, LH, and PDF contributed to acquisition of data, and all authors contributed to analysis and interpretation of data. EL performed statistical analysis. MER, LH, and PDF contributed to drafting the article. All authors revised the manuscript, provided intellectual input, and approved the submitted version.

Financial support. This study was funded by a contract from Magnolia Medical Technologies, Inc., to the University of Nebraska Medical Center.

The sponsor approved the final protocol but did not have a role in data acquisition, data analysis, or manuscript preparation.

Competing interests. MER has received past financial support from Magnolia Medical Technology, Inc., in the form of contracts between the sponsor and the University of Nebraska Medical Center for a prior prospective clinical trial and retrospective cohort analysis. In addition, MER has served as a consultant for Teleflex, 3M, and Citius Pharmaceuticals. All other authors report no potential conflicts of interest.

The authors thank Nicholas Conoan of the Electron Microscopy Core Facility (EMCF) at the University of Nebraska Medical Center for technical assistance. The EMCF is supported by state funds from the Nebraska Research Initiative (NRI) and the University of Nebraska Foundation and institutionally by the Office of the Vice Chancellor for Research.

References

1. Doern GV, Carroll KC, Diekema DJ, *et al.* Practical guidance for clinical microbiology laboratories: a comprehensive update on the problem of blood culture contamination and a discussion of methods for addressing the problem. *Clin Microbiol Rev* 2019;33:e00009–19. DOI: [10.1128/CMR.00009-19](https://doi.org/10.1128/CMR.00009-19).
2. Liaquat S, Baccaglini L, Haynatzki G, Medcalf SJ, Rupp ME. Clinical consequences of contaminated blood cultures in adult hospitalized patients at an institution utilizing a rapid blood-culture identification system. *Infect Control Hosp Epidemiol* 2021;42:978–984. DOI: [10.1017/ash.2022.22](https://doi.org/10.1017/ash.2022.22).
3. Alahmadi YM, Aldeyab MA, McElnay JC, *et al.* Clinical and economic impact of contaminated blood cultures within the hospital setting. *J Hosp Infect* 2011;77:233–236. DOI: [10.1016/j.jhin.2010.09.033](https://doi.org/10.1016/j.jhin.2010.09.033).
4. Ziegler MJ, Pellegrini DC, Safdar N. Attributable mortality of central line associated bloodstream infection: systematic review and meta-analysis. *Infection* 2015;43:29–36. DOI: [10.1007/s15010-014-0689-y](https://doi.org/10.1007/s15010-014-0689-y).
5. O'Grady NP, Alexander E, Alhazzani W, *et al.* Society of Critical Care Medicine and the Infectious Diseases Society of America guidelines for evaluating new fever in adult patients in the ICU. *Crit Care Med* 2023;51:1570–1586. DOI: [10.1097/CCM.0000000000006022](https://doi.org/10.1097/CCM.0000000000006022).
6. Sherertz RJ, Karchmer TB, Palavecino E, Bischoff W. Blood drawn through valved hub connectors carries a significant risk of contamination. *Eur J Clin Microbiol Infect Dis* 2011;30:1571–7. DOI: [10.1007/s10096-011-1262-6](https://doi.org/10.1007/s10096-011-1262-6).
7. Rupp ME, Cavalieri RJ, Marolf C, Lyden E. Reduction in blood culture contamination through use of initial specimen diversion device. *Clin Infect Dis* 2017;65:201–205. DOI: [10.1093/cid/cix304](https://doi.org/10.1093/cid/cix304).
8. Bates DW, Goldman L, Lee TH. Contaminant blood cultures and resource utilization. The true consequences of false-positive results. *JAMA* 1991; 265:365–369. DOI: [10.1001/jama.1991.03460030071031](https://doi.org/10.1001/jama.1991.03460030071031)
9. Souvenir D, Anderson DE Jr, Palpant S, *et al.* Blood cultures positive for coagulase-negative staphylococci: antiseptic, pseudobacteremia, and therapy of patients. *J Clin Microbiol* 1998;36:1923–1926. DOI: [10.1128/JCM.36.7.1923-1926.1998](https://doi.org/10.1128/JCM.36.7.1923-1926.1998).
10. Hankins R, Majorant OD, Rupp ME, *et al.* Microbial colonization of intravascular catheter connectors in hospitalized patients. *Am J Infect Control* 2019;47:1489–1492. DOI: [10.1016/j.ajic.2019.05.024](https://doi.org/10.1016/j.ajic.2019.05.024).